Titanized Mesh Implants

TiLENE® Plug Set

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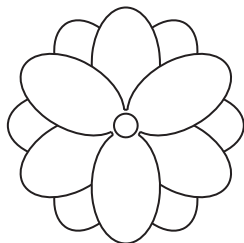
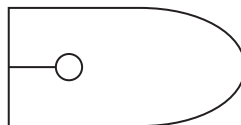
TiLENE® Onlay-Patch

Description

Large-pore, non-absorbable prosthetic mesh made from monofilament polypropylene thread with surrounding titanized coating. The mesh is perfectly coordinated with the dynamics of body tissue in terms of resistance to tearing and elasticity. The mesh implants are available in various basis weights with defined shapes and dimensions.

The TiLENE® Plug Set consists of one TiLENE® Plug and one TiLENE® Onlay-Patch.

These Instructions for Use are also valid for other TiLENE® Plug Set product variants.

Indication

Surgical mesh implants of titanized polypropylene mesh are used in the treatment of soft-tissue defects.

TiLENE® Plug Set is used to repair inguinal hernias.

Mode of action

The purpose of the implants is to strengthen tissue or bridge defects in soft tissue. Adsorption of a thin fibrous tissue layer, whose cells grow through the pores of the mesh fabric, results in the complete enclosure of the mesh material by the patient's own body tissue. Experimental investigations have shown that the implant remains soft and flexible in this context.

Warnings

- ☞ The following complications have been observed in the treatment of soft tissue defects with TiLENE® mesh implants: Seroma, hematoma, recurrence, infections, visceral adhesions, mesh migration and fistula formation.
- ☞ Although the mesh material is not affected by infections, the use of non-resorbable mesh material in an infected or contaminated wound always carries the risk of fistula formation and/or mesh extrusion.
- ☞ A foreign-body reaction triggered by the titanium coating may occur in the case of a titanium allergy.
- ☞ The titanium coating of the synthetic material will not affect radiological diagnostic techniques (X-ray, CT, MRI and ultrasonic examinations).
- ☞ Do not use the product if the packaging is open or damaged or should you have doubts about the sterility or if the use by date has expired.
- ☞ Each product is packed separately and is supplied EO-sterilised and pyrogen-free. It is only to be used once. Do not reuse, reprocess or resterilise. A reuse, reprocessing or resterilisation of disposable products can result in a reduced performance or a loss of their functionality. A reuse of disposable products can result in exposure to pathogens such as viruses, bacteria, fungi or prions.
- ☞ Damaged implants must not be used. Make sure that sharp objects do not affect the mesh beyond the usual clinical application during the operation.

Precautions

- In the case of wide dissection and large-scale tissue excision, it is recommended that hemostasis be carefully checked and any remaining cavity avoided wherever possible, so as to facilitate the rapid integration of hernia strengthening.
- This product should only be used by doctors who are familiar with surgical techniques of this type. The surgeon bears the responsibility for the implantation.
- The implant should only be handled with powder-free gloves to prevent the inflammatory effect of foreign particles such as talc.
- The implants are packaged in a sterile double-wrap packaging. Carefully check that the packaging is undamaged before use. Under no circumstances should the implant be used if the packaging is damaged, open or moist. The inner packaging should only be opened directly before insertion of the implant and should exclusively be handled with sterile gloves and sterile instruments.
- Before using TiLENE® mesh implants, the physician should check the size of the hernial orifice and ensure a sufficient overlap.

Recommended Implantation procedures

In the case of indirect hernias, the hernia sac is separated after preparation and ligature. To ensure optimal results select the right size of TiLENE® Plug. If necessary, trim inner ply of TiLENE® Plug in case less mesh material is required for smaller defects. Grasp the inner ply of TiLENE® Plug with a sterile forceps or clamp and push the plug gently through the internal ring, deep into the preperitoneal space. After positioning of the plug, secure the outer ply of the plug with 2-4 sutures.

Additionally TiLENE® Onlay-Patch is placed on the anterior surface of the posterior wall of the inguinal canal. The opening is placed around the spermatic cord to repair an additionally existing direct hernia and/or to provide strengthening to prevent the recurrence of a direct or indirect hernia.

Contraindications


All established standard contraindications and rules for soft-tissue implants apply to the use of TiLENE® mesh implants.

These include the following:

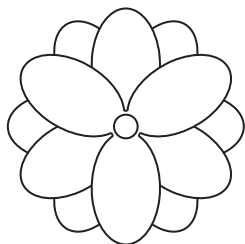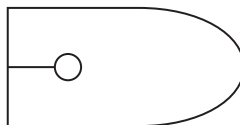
- Patients in the growth period (the known extensibility of this strengthening may be insufficient for the growth of children).
- Operation on a contaminated or inflamed area.

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